



SE-1/SE-3

Single/Three-channel ECG

Technical Specifications

	SE-1 Single Channel ECG	SE-3 Three-Channel ECG
Input Circuit	Defibrillator protection built-in	
Lead	Standard 12 leads	
Acquisition Mode	Single lead	12 leads simultaneously
A/D Converter	12 bits	
Measurement Range	± 5mVpp	
Time Constant	> 3.2s	
Frequency Response	0.05Hz ~ 150Hz	
Calibration Voltage	1mV ± 3%	
Sensitivity	2.5, 5, 10, 20 (mm/mV)	
Input Impedance	> 10MΩ	> 50MΩ
Input Circuit Current	< 50nA	
Noise Level	< 15 μVp-p	
Anti Baseline Drift	Auto	
Patient Leakage Current	< 10 μA (220V~240V)	
Patient Auxiliary Current	< 0.1 μA (DC)	
Filter	EMG Filter 35Hz (-3dB) AC Filter 50 Hz/60Hz (-20dB)	EMG Filter 25Hz/35Hz/45Hz /OFF DFT Filter 0.05/0.15/0.25/0.5Hz Lowpass Filter 150Hz/100Hz/75Hz AC Filter 50 Hz/60Hz
CMRR	>90dB	>110dB
Recording Mode	Auto: 1CH, 1CH+1R Manual: 1CH	Auto: 3CH, 3/2 CH, 3CH+1R, 1CH, 1CH+1R Manual: 3CH, 3/2 CH, 1CH Rhythm: Single rhythm lead
Rhythm Lead	Any lead selectable	
Safety Standard	IEC I/CF	
Power Supply	AC: 100V~115V/220V~240V, 50/60Hz Built-in rechargeable Li-ion battery: voltage=14.4V More than 3 hrs continuous work under battery supply	
Recorder	Thermal printer	
Record Paper	50mm x 30m	80mmx20m
Paper Speed	25mm/s, 50mm/s (± 3%)	5mm/s, 6.25mm/s, 10mm/s, 12.5mm/s, 25mm/s 50mm/s (± 3%)
Dimension	288mm x 210mm x 70mm	
Weight	2.3kg (with battery)	2.5kg (with battery)
Display	192 x 64 pixels LCD	320 x 240 pixels LCD (Wide Screen) 192 x 64 pixels LCD (Narrow Screen)
External Input/Output	Input: > 100 kΩ ; Sensitivity 10mm/V ± 5% Output: < 100 Ω; Sensitivity 1V/mV ± 5%	
Communication Interface	RS232	Ethernet, RS232

Features and Specifications are subject to change without prior notice



SE-1/SE-3

Single/Three-channel ECG



Edan Instruments, Inc.

3F-B, Nanshan Medical Equipments Park, Nanhai Rd 10199, Shekou, Nanshan, Shenzhen, 518067 P.R. China
Tel: 86-755-2689 8326, 2689 8316 Fax: 86-755-2688 2223 www.edan.com.cn Email: info@edan.com.cn



KABIOS suarl - Rue Hateb Ibn Abi Baltâa 4023 Erriadh - Sousse - Tunisie

Tél: (216) 73 30 70 44 - Fax: (216) 73 30 70 45

E-mail: contact@kabios-biomedical.com - Web: www.kabios-biomedical.com



SE-3 Three - Channel ECG



SE-3(Wide Screen)

- Foldable LCD screen, vivid waveform display
- 12 leads simultaneous acquisition
- Built-in high-resolution thermal printer
- Built-in rechargeable Li-ion battery, AC/DC power supply
- Automatic ECG measurement and interpretation

Smart ECG Viewer

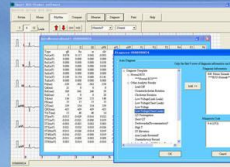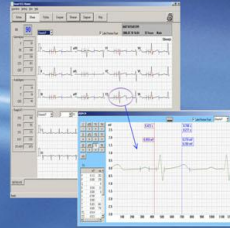



PC-based ECG data management software

- Complete digital filtering, resisting baseline drift, AC and EMG interference
- Automatic baseline adjustment
- Auto/manual work mode
- Lead-off detection and alert
- Three printing formats (Auto, manual and rhythm) selectable
- Heart rate variability (HRV) analysis



SE-3(Narrow Screen)

- 12 leads printing by external jet-printer
- Large patient data storage: 120 ECGs (Narrow Screen) 144 ECGs (Wide Screen)
- PC-based ECG data management software (option)
- ECG data input function support

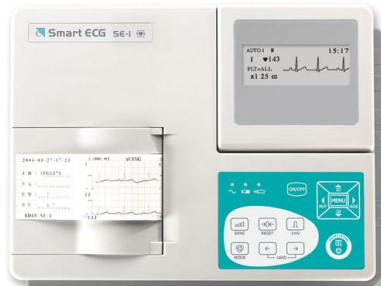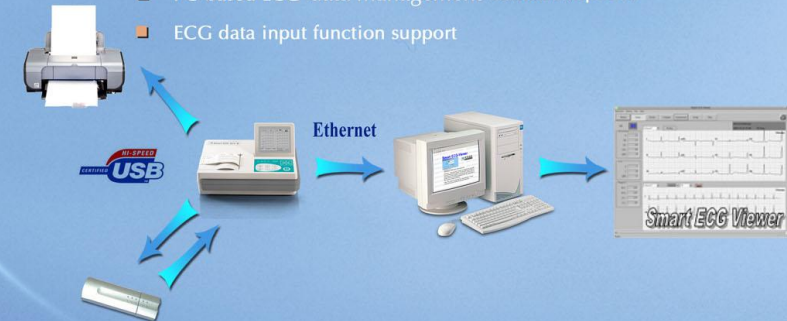

SE-1 Single Channel ECG

- Foldable LCD screen, vivid waveform display
- Built-in high-resolution thermal printer
- Built-in rechargeable Li-ion battery, AC/DC power supply
- Well designed membrane keyboard
- Digital filtering, resisting baseline drift, AC and EMG interference
- Automatic baseline adjustment
- Auto/manual work mode
- Lead-off detection and alert

